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Read this manual carefully before use, and retain it for future reference.

Safety Instructions






When using the toothbrush for the first time, it may cause itch or minor gum bleeds. This is normal when the gums are not accustomed to the efficient cleaning of an electric toothbrush. It is recommended to use gentle mode when your gums are not used to it yet. Normally, gum bleeds will stop within one to two weeks. Stop using the toothbrush and consult a dentist if it persists.

⚠ WARNINGS

- The plug must be removed from the socket-outlet before cleaning or maintaining the appliance.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the appliance.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Do not expose the appliance or battery to excessive temperatures.
- Be aware of the risk of terminals of the battery-operated appliance or battery being short-circuited by metal objects.
- This appliance contains a battery that is non-replaceable.
- Only use and charge the appliance at temperature between 5°C and 35°C.

Note: Illustrations of product, accessories, and user interface in the user manual are for reference purposes only. Actual product and functions may vary due to product enhancements.

Product Overview

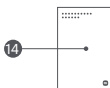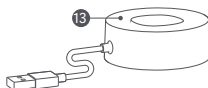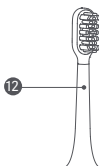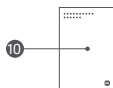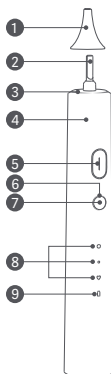- | | |
|-----------------------------------|----------------------------------|
| 1 Shaft Protective Cover | 9 Battery Level Indicator |
| 2 Shaft | 10 User Manual |
| 3 Colorful Smart Indicator | 11 Toothbrush Head Cover |
| 4 Toothbrush Handle | 12 Toothbrush Head |
| 5 Power Button | 13 Charging Dock |
| 6 AI Ring Indicator | 14 Warranty Notice |
| 7 Mode Button | |
| 8 Mode Indicators | |
| ○ Standard Mode | |
| ○ Gentle Mode | |
| ♥ Custom Mode | |


How to Use

1. Installing the toothbrush head


- Make sure that the toothbrush head is installed in the same direction as the front of the toothbrush handle, and make sure there is a slight gap in between the toothbrush head and the handle. The gap is required to avoid collision of the toothbrush head with the handle during high-frequency oscillation, so as to generate more power from the motor.

2. Wet the toothbrush head and squeeze a proper amount of toothpaste onto the toothbrush head.

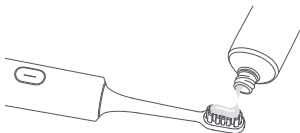3. Turning on/turning off/switching modes

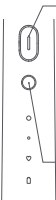Turn on: Press the power button to turn the toothbrush on.

Turn off/pause: Within ten seconds of the toothbrush being turned on, press the power button to turn it off. Ten seconds after the toothbrush is turned on, press the power button to pause it. If there is no operation within 30 seconds of being paused, the toothbrush will automatically turn off.

Switch modes: Press the mode button to cycle through different modes. Once switched, the corresponding mode indicator will light up.

Bluetooth broadcast: When the toothbrush is turned off, press and hold the mode button for over seven seconds until the colorful smart indicator breathes orange, indicating that the Bluetooth broadcast is enabled; press and hold the mode button again for over seven seconds until the colorful smart indicator keeps orange briefly and then turns off, indicating that the Bluetooth broadcast is disabled.

Restore: Press and hold the power button and mode button at the same time for over seven seconds to restore the toothbrush to factory settings, disconnect it from the app, and clear all app settings.

Reset: Press and hold the power button for over 12 seconds to restore the toothbrush to factory settings and clear all app settings, but keep it connected to the app.

4. Healthy Brushing

- Hold the toothbrush head at an angle of 45° to the gum line, and slowly move it along the teeth and gum line.
- While brushing, the toothbrush pauses briefly every 30 seconds, to prompt the user to change the brushing position.
- When there is no active pause during use, the toothbrush will automatically turn off once the brushing mode ends.

- Three brushing modes are available:

Standard mode: Suitable for users with good oral health.


Gentle mode: Recommended for users with oral sensitivity.

Custom mode: A two-minute beginner mode provided by default. More brushing modes are available after connecting the toothbrush to the Mi Home/Xiaomi Home app.

Note: The supplied toothbrush heads are whitening toothbrush heads with strong cleaning power. Therefore, it is recommended that users with any oral sensitivity should use the gentle mode.

5. Cleaning the toothbrush head

- After brushing, rinse the toothbrush head with running water.
- After rinsing, put on the toothbrush head with the protective cover to keep it clean and hygienic.

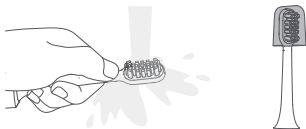Note: It is recommended to replace the toothbrush head every three months.

Smart Functions

- The toothbrush is equipped with an overpressure protection system that detects excess pressure in real time. If too much pressure is applied when you brush your teeth, the colorful smart indicator will keep red and the toothbrush will vibrate less to remind you.

- Connect the toothbrush to the Mi Home/Xiaomi Home app to set the brushing mode, enable the smart AI adjustment, view your brushing data and oral health report, enable the spill-resistant feature, set custom colors for the colorful smart indicator, and enjoy other personalized features.
- When the toothbrush has power, picking it up will wake it up. After connecting it to the app and enabling the shake function, shaking the toothbrush quickly twice when it is awake will change the color of the colorful smart indicator.

CAUTION: When the toothbrush is turned on, it cannot be connected to the app.

Connecting with the Mi Home/Xiaomi Home App

- This product works with the Mi Home/Xiaomi Home app*. Control your device with Mi Home/Xiaomi Home app.
- Scan the QR code to download and install the app. You will be directed to the connection setup page if the app is installed already. Or search for "Mi Home/Xiaomi Home" in the app store to download and install it.

- Open the Mi Home/Xiaomi Home app, tap "+" on the upper right, and then follow the instructions to add your device.

9FAB2DEF

- * Note: The toothbrush's Bluetooth broadcast is disabled by default and should be enabled before the toothbrush is connected to the Mi Home/Xiaomi Home app. When the toothbrush is turned off, press and hold the mode button for over 7 seconds until the colorful smart indicator breathes orange, indicating that the Bluetooth broadcast is enabled.
- * The app is referred to as Xiaomi Home app in Europe (except for Russia). The name of the app displayed on your device should be taken as the default.

Note:

The version of the app might have been updated, please follow the instructions based on the current app version.

- Colorful smart indicator status description

Toothbrush status	Colorful smart indicator status
Pick-up-to-wake	Breathing for nine seconds before turning off (the default color is sky blue, which can be customized in the app)
Paused	Breathing blue
Overpressure alert	Keeping red

Colorful Smart Indicator

AI Ring Indicator

- AI ring indicator status description


Toothbrush status	AI ring indicator status
Smart AI adjustment disabled	When the toothbrush is turned on and in use, the AI ring indicator keeps off.
Smart AI adjustment enabled	When the toothbrush is turned on and in use, the AI ring indicator keeps blue.

Charging

- Battery level indicator status description

Toothbrush status	Battery level indicator status	Battery level
Using	Off	20%–100%
	Breathing orange	10%–20%
	Blinking orange three times (the toothbrush cannot be turned on)	<10%
Charging	Off	100%
	Breathing white	20%–100%
	Breathing orange	<20%

- The charging dock of the toothbrush is designed for the T501 and is not compatible with other products.
- When the battery level indicator blinks orange three times and the toothbrush won't turn on, it will take about 6.5 hours to fully charge the toothbrush (the exact time depends on the amount of battery left).

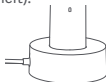- Please use a 5 V $\overline{\text{---}}$ 1 A charger (sold separately) when charging.

Troubleshooting

Issue	Possible Causes	Solutions
Unable to turn on by pressing button	Charging	Remove the toothbrush from the charging dock
	Battery level is too low	Place the toothbrush on the charging dock to charge
Unable to install the toothbrush head onto the shaft	There is something stuck inside the toothbrush head	Use a tool to remove the object from the toothbrush head
	There is something stuck on the shaft	Clean the shaft
	Not an original toothbrush head	Replace it with an original toothbrush head
Excessive power consumption and charging frequently	Paired to the phone via Bluetooth for a long time	Disable Bluetooth when not using the toothbrush
	Used too frequently	Place the toothbrush on the charging dock to charge
	The issue is still not resolved after excluding the causes above.	Contact the after-sales service team.
Error with the vibration of the toothbrush head	Vibration varies depending on the modes set in the app.	This is normal. Switch modes as desired.
	toothbrush head is not installed properly	Reinstall the toothbrush head
	There is something stuck between the toothbrush head and the handle.	Remove the toothbrush head and remove the object
Unable to connect to the app	The Bluetooth of your phone is not enabled	Enable Bluetooth on your phone, and then follow the instructions in the app to connect.
	The toothbrush's Bluetooth broadcast is disabled.	When the toothbrush is turned off, press and hold the mode button for over seven seconds until the colorful smart indicator breathes orange to enable the Bluetooth broadcast. Then try to reconnect the toothbrush to the app.

Unable to connect to the app	Toothbrush is too far away from the phone	Place the phone near the toothbrush
	Toothbrush is vibrating	Turn off the toothbrush and then connect again.
	App is lagging.	Close the app and then reconnect.
	The toothbrush crashes.	Press and hold the power button for over 12 seconds to reset the toothbrush, and then reconnect.

Specifications

Name: Electric Toothbrush	
Model: MES607	
Rated Voltage: 3.6 V ---	Rated Power: 2.5 W
IP Rating: IPX8	Battery: 3.6 V --- , 2550 mAh, 9.18 Wh
Charging Time: Approx. 6.5 hours	
Compatible with: Android 8.0 & iOS 12.0 or above	
Wireless Connectivity: Bluetooth® 5.0 BLE (Bluetooth Low Energy)	
Maximum Output Power: < 20 dBm	
Operating Frequency: 2400–2483.5 MHz	

Name: Charging Dock	
Model: MES607	Rated Input: 5 V --- 1 A
Charging Method: Wireless charging	
Maximum Output Power: -20 dBm	
Operating Frequency: 110–205 KHz	
This device has been tested and meets applicable limits for radio frequency (RF) exposure.	

Environmental Protection


When the toothbrush reaches the end of its service life, please remove, dispose of and recycle the built-in battery according to local laws and regulations of the country or region where the product is used. If the battery cannot be removed, dispose of it in strict accordance with the regulations on classified domestic waste disposal in your area.

WARNINGS


- Do not expose the battery to fire to prevent the battery from exploding or releasing toxic and hazardous substances.
- Do not short-circuit the battery to prevent it from combusting and causing a fire.
- Do not disassemble the toothbrush unless you want to dispose of it; otherwise, the toothbrush will be damaged.
- When the battery is at end of life, the appliance shall be properly disposed of.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.

Removing Battery


- Insert a flathead screwdriver into the gap around the end of the handle to pry off the bottom cover.

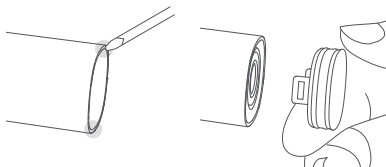- Use a flathead screwdriver to press the two clips symmetrically arranged on both sides while pressing the shaft at the top of the toothbrush against a hard flat surface, to release the inner compartment.

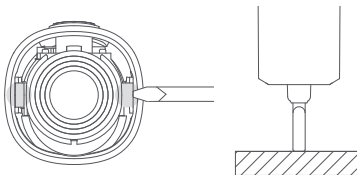- Cut the three wires connecting the battery with the PCB board.
- Pry off the battery from the back of the PCB board, and dispose of it at an officially designated recycling center.

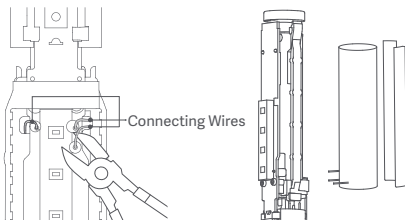WEEE Disposal and Recycling Information

All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste. Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by the government or local authorities. Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact the installer or local authorities for more information about the location as well as terms and conditions of such collection points.

Manufactured for: Xiaomi Communications Co., Ltd.

Manufactured by: Shenzhen Soocas Technology Co., Ltd.
(a Mi Ecosystem company)

Address: Room 1101, 1102, 1103, 1104, Building 2 Chongwen Area, Nanshan iPark, No. 3370 Liuxian Avenue, Fuguang Community, Taoyuan Subdistrict, Nanshan District, Shenzhen, China

For further information, please go to www.mi.com

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WARRANTY NOTICE

As a Xiaomi consumer, you benefit under certain conditions from additional warranties. Xiaomi offers specific consumer warranty benefits which are in addition to, and not instead of, any legal warranties provided by your national consumer law. The duration and conditions related to the legal warranties are provided by respective local laws. For more information about the consumer warranty benefits, please refer to Xiaomi's official website <https://www.mi.com/global/support/warranty/>. Except as prohibited by laws or otherwise promised by Xiaomi, the after-sales services shall be limited to the country or region of the original purchase. Under the consumer warranty, to the fullest extent permitted by law, Xiaomi will, at its discretion, repair, replace or refund your product. Normal wear and tear, force majeure, abuse or damage caused by the user's negligence or fault are not warranted. The contact person for the after-sale service may be any person in Xiaomi's authorized service network, Xiaomi's authorized distributors or the final vendor who sold the products to you. If in doubt please contact the relevant person as Xiaomi may identify. The present warranties do not apply in Hong Kong and Taiwan.

Products which were not duly imported and/or were not duly manufactured by Xiaomi and/or were not duly acquired from Xiaomi or a Xiaomi's official seller are not covered by the present warranties. As per applicable law you may benefit from warranties from the non-official retailer who sold the product. Therefore, Xiaomi invites you to contact the retailer from whom you purchased the product.

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